




MIPAR™

Release Notes

2.1.1 - 2.1.4

Color Select

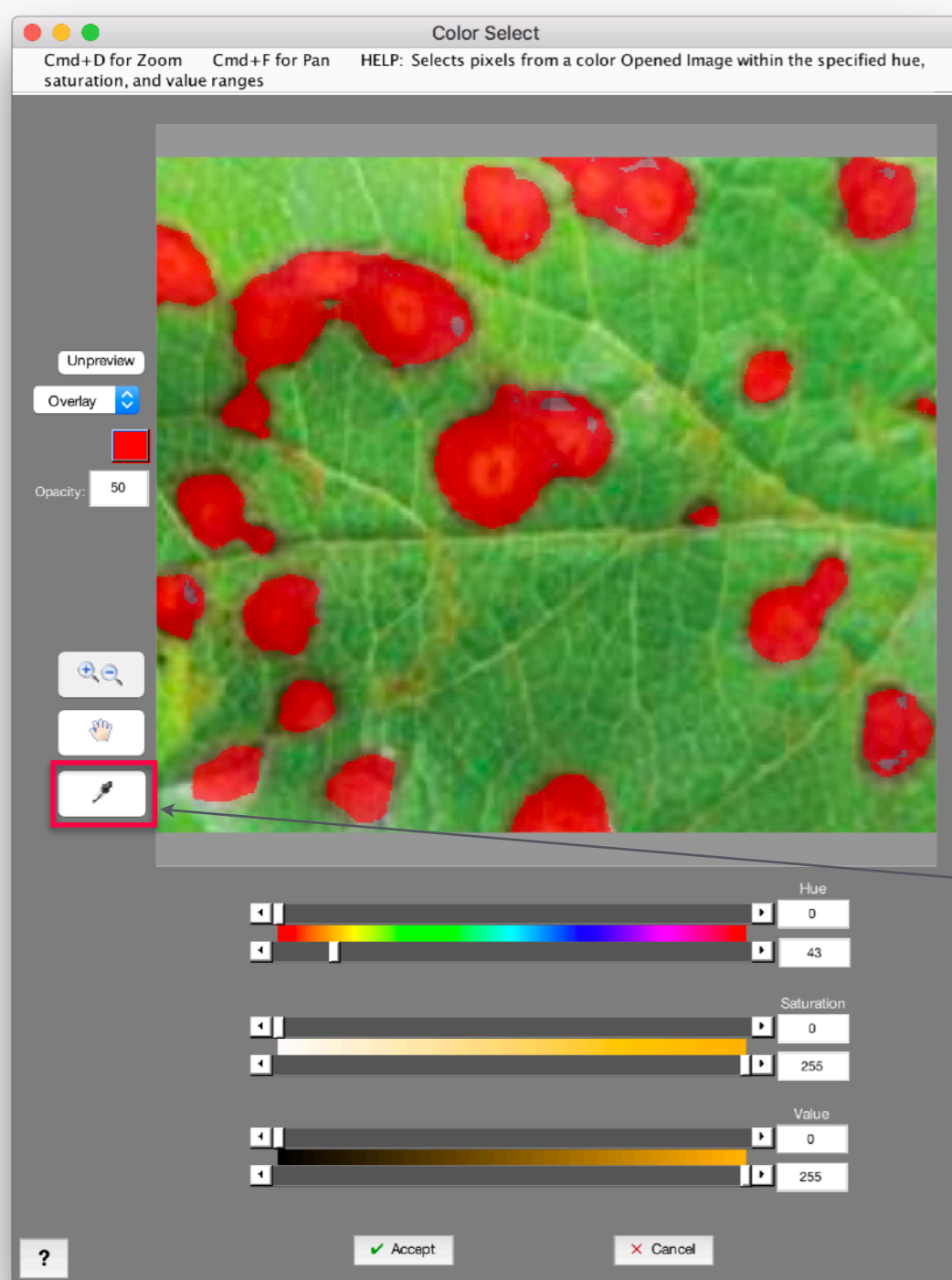

Easy Color Selection

Drag HSV sliders to zero-in on the precise color of interest

Point-and-Click

With the "Choose Color" tool, click a spot in the image, and the HSV limits are set for you

Bug Fixes and Other Changes

- Several fixes and improvements
- *Manual Edit* steps now remember your last settings
- New Edit > *Draw Random Lines* lets you place lines for visual intercept measurements
- *Color Cluster*, *Color Select*, and *Channel Operation* are now different step classes
- *Watershed* window now displays previous image as the reference rather than original
- Fixed issue with layer name printed during batch feature measurements in Post Processor